

vygonvet

PRODUCT FOCUS



WHY CHOOSE VYGON?



**OVER
60
YEARS
EXPERIENCE**

PATIENT CARE LIES AT THE HEART OF WHAT WE DO



**CLINICAL
EXPERTISE**



**PLANET
MARK
CERTIFIED
BUSINESS**

**MORE THAN
100 MILLION
PRODUCTS IN
110 COUNTRIES**



CONTENTS

Our commitment to the environment	3
Needle-free devices	4
IV management	8
CVCs	10
Regional nerve blocks	12
Cannulas	14
MostCare Up	15
Custom packs	16
Oncology	18
Oncology PPE	19
Nutrisafe2® - safe enteral feeding system	20
Enteral feeding	21
Theatre	24

OUR SERVICE OFFERING



Customer Service & Technical Support

Talk through your enquiry with our dedicated Teams.



Education & Training

We offer our customers a variety of valuable and comprehensive training options to help you and your teams meet your training requirements.

This Product Focus features a selection of our products.

Visit **vygonvet.co.uk** for our full range.



A leading supplier of single-use veterinary consumables for more than 20 years, Vygon Vet has built an enviable reputation for high quality products and excellent customer service.

Working closely with practice teams focused on delivering the best care for their patients, our priority is to tailor our products and support to ensure surgical procedures run smoothly and efficiently, treatments and therapies deliver faster with better recovery. This philosophy is behind the success of all our products including the popular needle-free IV access system Bionector and Custom packs.

In addition, we also focus on supporting practices to push through treatment boundaries and benefit from innovation. Ground-breaking products include the Plexygon Nerve Stimulator, to access the patient benefits of regional anaesthesia, and Qimono, our closed system transfer device, which provides protection for veterinary teams whilst delivering life-saving chemotherapy.

Vygon Vet is a division of leading healthcare and medical products supplier Vygon (UK) Ltd.

OUR COMMITMENT TO THE ENVIRONMENT

2021 was a landmark year as Vygon UK achieved carbon neutrality in accordance with the guidance set out in PAS 2060, with certification renewed in 2022 as year-on-year emission reductions were achieved and the residual emissions were offset with the purchase of high-quality Verified Carbon Standard (VCS) emission reduction projects. Vygon are in the process of expanding its emissions reporting to include Scope 3 emissions, this will involve a baseline year reset and result in reporting higher emissions than in previous years and therefore, unable to purchase carbon credits. Following the baseline year and emission targets reset, our focus will be to regain Carbon Neutral status.

We are delighted to announce that Vygon UK has successfully undergone the NHS Evergreen Sustainable Supplier Assessment, attaining level 2 status. This is key to supporting us to understand how we align to the NHS's long-term sustainability priorities and the pathway to progress.

We're proud to be a Planet Mark Certified Business, marking a vital first step on our journey to net zero. We annually measured our carbon emissions with Planet Mark, as we cannot manage what we do not measure. Vygon UK are committed to a 5% annual reduction in Scope 1 and 2 emissions. We look forward to sharing our progress through Planet Mark's Net Zero Certification Programme.

These accomplishments underscore our unwavering dedication to sustainable practices and reinforces our role as a responsible contributor to the healthcare system.



Discover our journey
scan with your smart device




NEEDLE-FREE DEVICES

bionector®

Bionector delivers safe, needle-free IV access for sampling, injections, and infusions. Its self-sealing, easy-to-disinfect membrane and closed design protect both patients and clinicians. Fully Luer-lock compatible, Bionector offers reliable performance across your entire infusion workflow.


Code	Description	Box qty
00089603	bionector closed system needle-free device	50

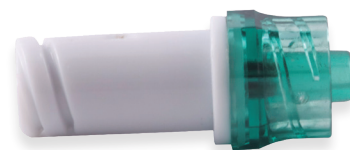

bionector® T-Piece

Allows intermittent injection via latex-free membrane during simultaneous infusion.

Code	Description	Box qty
1VTUSBN-448	10cm luer-lock T-piece extension connector with bionector	50

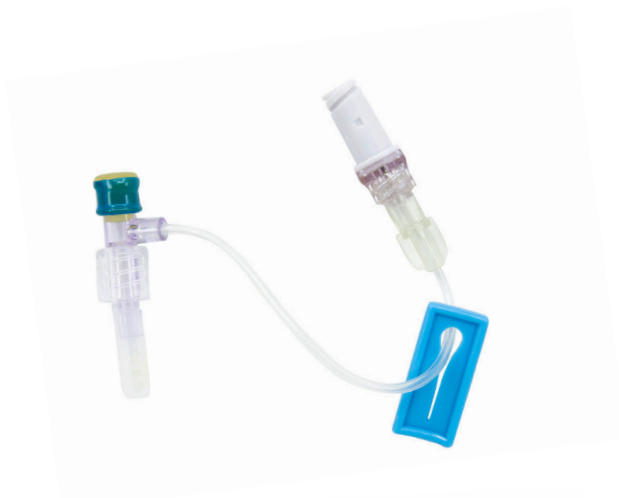

bionector® T-piece with two needle-free access ports

Totally needle-free T-piece device, reduces risk of needle-stick injuries to both patients and staff, as well as reducing the risk of infection.

Code	Box qty
1VTUSBN-450	50



bionector® vial caps

For use when aspirating drugs from capped vial and drug reconstitution.


Code	Size (mm)	Box qty
00082101	14	50
00082102	20	50



bionector® and bag spike

Code	Box Qty
00085611	50



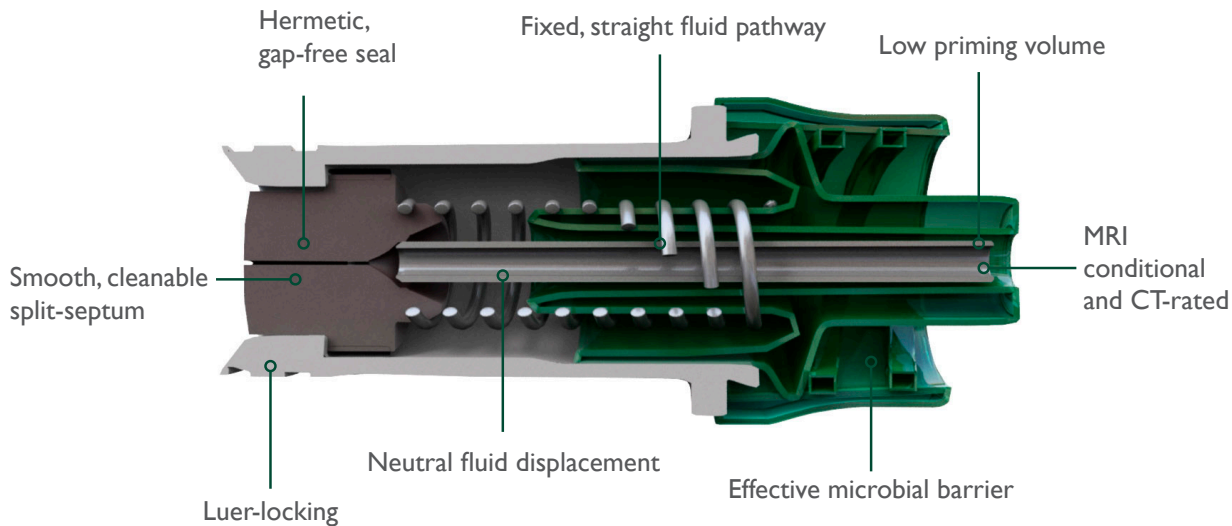

Exceeding global opinion leaders' recommendations

An established market-leading needle-free device which meets the full range of global opinion leaders' recommendations for reducing CRBSIs. It has been proven to provide an effective barrier against microbial ingress and help standardise practice by combining a fixed, straight fluid pathway with innovative neutral displacement technology.

MRI conditional and CT-rated

Bionector is proven not to represent any risk to either patients or practitioners during an MRI of up to three Teslas. CT-rated for use with power injectors, Bionector has a maximum pressure resistance of 350psi and a maximum flow rate of 10mL/s.

Neutral displacement

Leading the way with a neutral fluid displacement, a specific post-flushing clamping sequence is not required. This helps prevent blood reflux and reduce catheter occlusions.

Clinically proven

Backed up by a robust library of clinical studies Bionector is proven to be easy to clean and clear. Its smooth split-septum fits tightly into the device housing ensuring it is free from any gaps. The straight, fixed fluid pathway has been proven 'easy to clear', designed to provide the most direct and least tortuous route with no moving parts (such as mechanical valves), which reduces the surface area available for biofilm formation.

Low deadspace

The straight fluid pathway is proven 'flushable' for macro and microscopic particles such as blood. This is due to a minimal deadspace of just 0.018mL allowing for a low flushing volume (5mL) to clear the device.



bionector® octopus double lumen

Octopus needle-free extension set, including double lumen Bionector octopus extension with 10cm tube and clamp PUR.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
000841264	10	100	1.5	2.5	50



bionector® octopus double lumen with two anti-reflux valves

Octopus needle-free extension set, including double lumen Bionector octopus extension and two anti-reflux valves with 0.4mL capacity.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
000841232	3	120	1.5	2.5	10



bionector® octopus triple lumen

Octopus needle-free extension set, including triple lumen Bionector octopus extension with 10cm tube and clamp PUR.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
000841364	10	100	1.5	2.5	10



bionector® octopus triple lumen with three anti-reflux valves

Octopus needle-free extension set, including triple lumen Bionector octopus extension and three anti-reflux valves with 0.4mL capacity.

Code	Length (cm)	Flow rate (mL/min)	ID (mm)	OD (mm)	Box qty
000842312	3	120	1.5	2.5	10



3M™ Curot™ disinfecting port protectors

A disposable passive disinfection device. Curot port protectors ensures needle-free access devices are always clean, therefore removing the need to actively disinfect before and after each access.

Code	Description	Box qty
VWJIM001-270RB	Curot™ disinfecting port protector compact patient box	30

Curot™ is a trademark of 3M Company.



3M

Authorised Distributor
Medical Products




IV MANAGEMENT

lectrocath

lectrocath 50cm micro bore male/female luer-lock PE extension set.

Lectrocath delivers secure, chemically inert polyethylene extension lines for IV or arterial catheters, offering Luer-lock compatibility, high pressure resistance, and easy colour-coded identification.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
00115505	0.7	2.0	50	1.0	50

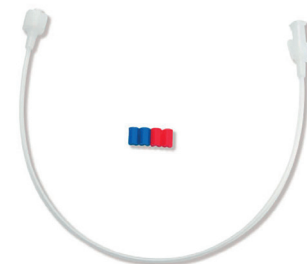

lectrospiral

Coiled IV extension line in polyurethane with luer-lock connectors.


Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
00115560	100	1.2	1.0	2.0	15
00115565	150	1.6	1.0	2.0	15
00115570	200	2.0	1.0	2.0	15
00115580	300	2.9	1.0	2.0	15
00115590	400	3.8	1.0	2.0	15
00115960	100	2.0	1.5	2.5	15
00115965	150	2.8	1.5	2.5	15
00115970	200	3.6	1.5	2.5	15
00115980	300	5.2	1.5	2.5	15
00115990	400	6.9	1.5	2.5	15



v-green IV extension line

PE extension line with male to female luer-lock connections. Narrow bore to minimise priming volume, ideal for drug infusion in animals of all sizes.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
07110005	50	0.7	1.0	2.0	50
07110010	100	1.2	1.0	2.0	50
07110015	150	1.6	1.0	2.0	50
07110020	200	2.0	1.0	2.0	50
07110105	50	1.1	1.5	2.5	50
07110110	100	2.0	1.5	2.5	50
07110115	150	2.8	1.5	2.5	50
07110120	200	3.6	1.5	2.5	50

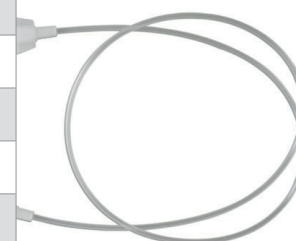

vyctic three-way stopcock

Lipid resistant, lockable three-way stopcock in PC-PE, with one rotating male luer and two female luer-lock ports. Colour-coded plugs for venous or arterial identification, and the handle position indicates open ports.

Code	Description	Box qty
00087620	Lockable three-way stopcock in PC-PE	50



protect-a-line extension line

Triple-layer extension line with anti-siphon valve and clamp—engineered to prevent unintended flow and ensure controlled, secure syringe-pump infusion.

Code	Length (cm)	Capacity (mL)	ID (mm)	OD (mm)	Box qty
00083201	150	1.4	1.0	2.5	50
000832211	200	1.3	1.5	3.0	50



CVCS

multicath® forward

Radiopaque multi-lumen central venous catheter with flexible tip. 1cm graduations start at 9cm from the distal tip. Each lumen has a different hub colour to aid identification. The guidewire feeder aids insertion whilst still allowing the user to sense any resistance. Each kit includes an introducer needle and a short PTFE cannula to give the user clinical choice for insertion.

Contents:

- 1 x multicath forward catheter
- 1 x Introducer needle 18G x 70mm
- 1 x Mosquito introducer cannula 18G x 64mm
- 1 x "J" guidewire in feeder, length 60cm, OD 0.9mm
- 1 x Vessel dilator
- 2 x Latex-free injection membranes
- 1 x Mini scalpel N
- 1 x 5mL luer-slip syringe
- 1 x Secondary fixation wing and anchor

Miniaturised colour-coded hubs

- Minimise dead space
- Maximise comfort
- Allows easy identification of infusions

Extension lines

- To manipulate away from the exit site

Fixation wings

- Small and soft for fragile skin

Choice of catheter lengths and sizes:


Offering practitioners, a choice when selecting a catheter to suit the patient and therapy needs.




3Fr




4.5Fr

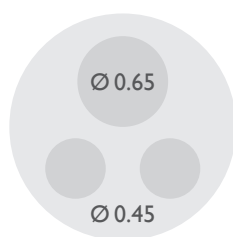

4.5Fr

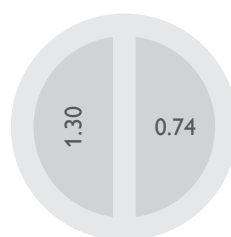

5.5Fr

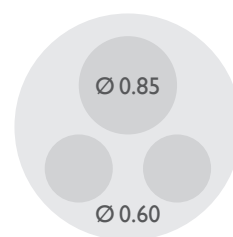

5.5Fr

multicath® forward

Code	Lumen	Catheter (PUR)			Flow rate & priming volume (mL/min)		
		Length (cm)	Size (Fr)	OD (mm)	Dist	Med	Prox
001202063	2	6	4.5	1.5	22 0.20	-	24 0.20
001203082	2	8	4.5	1.5	19 0.22	-	22 0.22
001203123	2	12.5	4.5	1.5	15 0.25	-	18 0.25
000157508	2	8	5.5	1.85	46 0.34	-	77 0.38
000155508	3	8	5.5	1.8	27 0.22	16 0.21	17 0.24
000155512	3	12.5	5.5	1.8	24 0.24	12 0.23	13 0.23

smartmidline

smartmidline is our CT-rated midline (up to 325psi and 7mL/s). It is simple and quick to place using the direct Seldinger technique. Suitable for a wide range of clinical indications from 2Fr to 5Fr and 4cm to 25cm, it is also an integrated extension line, which helps to reduce the risk of post insertion complications.

Code	Size (Fr)	Length (cm)	Guidewire (cm)	Max flow rate* (mL/s)	Box qty
0001281208	2	8	26	1.0	20
0001281310	3	10	26	1.5	20
0001281315	3	15	40	1.0	20
0001281415	4	15	50	5.0	20

*Viscosity 11.8cP



CVC Webinar

European Specialist in Veterinary Anaesthesia and Analgesia pinpoints potential mistakes and varying techniques of CVC placement in this video guide. Scan the QR code for more details.




REGIONAL NERVE BLOCKS

echoplex needles

Designed for clinicians who use ultrasound and stimulation techniques to locate the peripheral nerve sheath.

Code	Size (G)	Length (mm)	Box qty
006194103	21	100	20
006194123	21	120	20
006194153	21	150	20
006194253	23	25	20
006194353	23	35	20
006194503	22	50	20
006194853	21	85	20

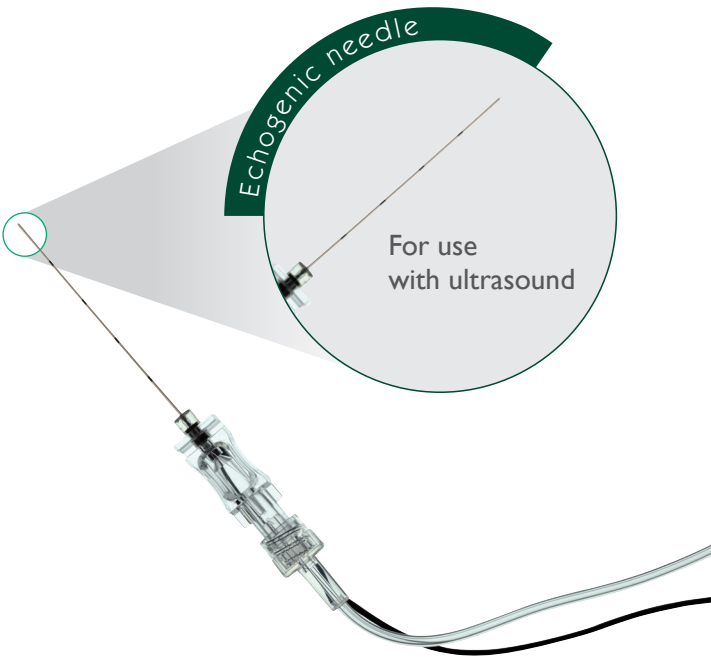

plexygon

Nerve stimulator allowing nerve location through the emission of a low intensity electrical current.

Code	Description	Box qty
00750131	plexygon nerve stimulator	1
08750302A	Replacement cable (Red and black)	1



Plexygon Nerve Stimulator Webinar

European Specialist in Veterinary Anaesthesia and Analgesia pinpoints the benefits of using the Plexygon Nerve Stimulator in this video guide. Scan the QR code for more details.



“We have experienced fewer problems with blood pressure, and we require lower levels of opiates. The dog can reliably go home the same day because the effects will last for around 8 to 12 hours.”

Mark Holmes
Clinical Director, Avondale Vets

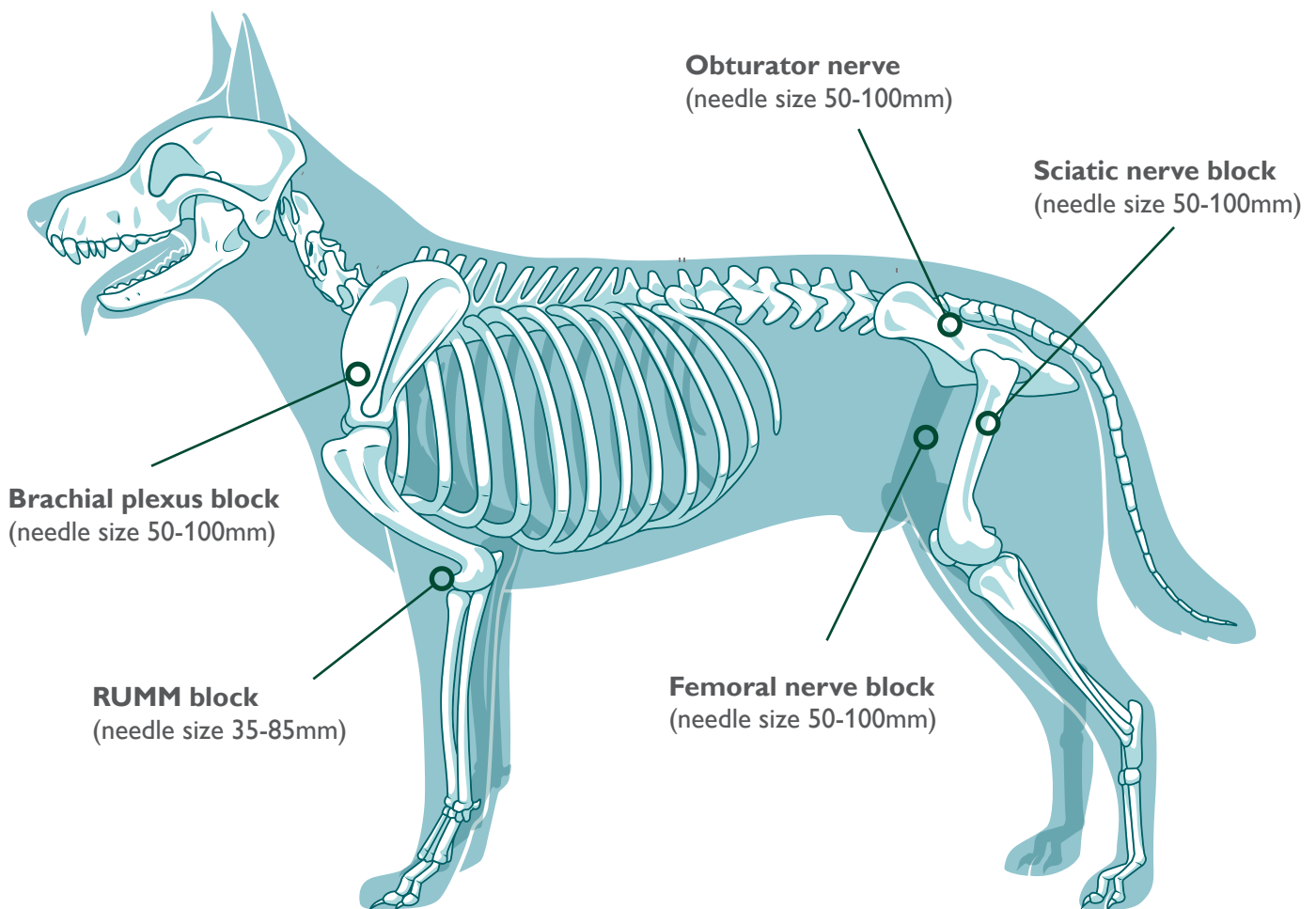


Image shown is for illustration purposes only.






**Speak to a product specialist about trialling a device.
Various purchase options are available.**

CANNULAS




intraflon non-ported

A short IV cannula in PTFE with transparent ultra-thin walls to optimise flow rate.

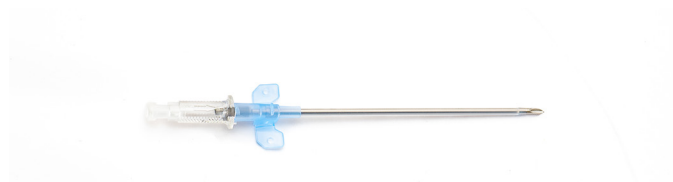

Code	Colour	Size (G)	ID (mm)	OD (mm)	Length (cm)	Max flow rate* (mL/s)	Box qty
00012106		24	0.4	0.7	18	18	50
00012108		22	0.5	0.8	28	24	50
00012110		20	0.7	1.0	30	52	50
00012113		18	0.9	1.3	38	100	50
00012116		16	1.2	1.6	52	174	50
00012121		14	1.6	2.1	52	260	50
00012127		12	2.2	2.7	52	450	50

intraflon non-ported longer

Code	Colour	Size (G)	ID (mm)	OD (mm)	Length (cm)	Max flow rate* (mL/s)	Box qty
00012210		20	0.7	1.0	36	50	50
00012213		18	0.9	1.3	47	95	50
00012216		16	1.2	1.6	60	170	50
00012221		14	1.6	2.1	80	250	50
00012227		12	2.2	2.7	80	430	50



Haemodynamic Monitoring

MostCare Up is haemodynamic monitoring at its most effective. With accurate patient monitoring minute by minute, beat to beat. With its intuitive, simple to use and customisable interface, you can set up and monitor a patient's cardiac output in 60 seconds, with no need for calibration.

MostCare Up uses a robust patented algorithm known as the Pressure Recording Analytical Method (PRAM™), giving real time data to help you make the right decisions at the right time. Our endless system has unlimited patient usage and does not require consumables, ensuring you can optimise its usage without additional cost and helps improve sustainability. MostCare Up can also connect to your patient information systems to help you stay more efficient.

Research article 'Briganti et al A preliminary study evaluating cardiac output measurement using Pressure Recording Analytical Method (PRAM) in anaesthetized dogs' available upon request.



Fast

60 second setup



Simple

No calibration required
Intuitive, customisable interface



Accurate

Robust PRAM™ algorithm
Patented dynamic filter to enhance signal quality



Lifetime value

Unlimited patient usage
No consumables required,
improving sustainability



Data management

HL7 compatible to enable data transfer




**Speak to us about trialling a device.
Various purchase options are available.**

CUSTOM PACKS

A custom pack allows you to create a high quality pack to match your local protocol, incorporating all the items you need.

Custom packs promote standardised insertion procedures, ensuring maximum barrier precautions and reduced risk of infection, ideally suited for many therapy areas.

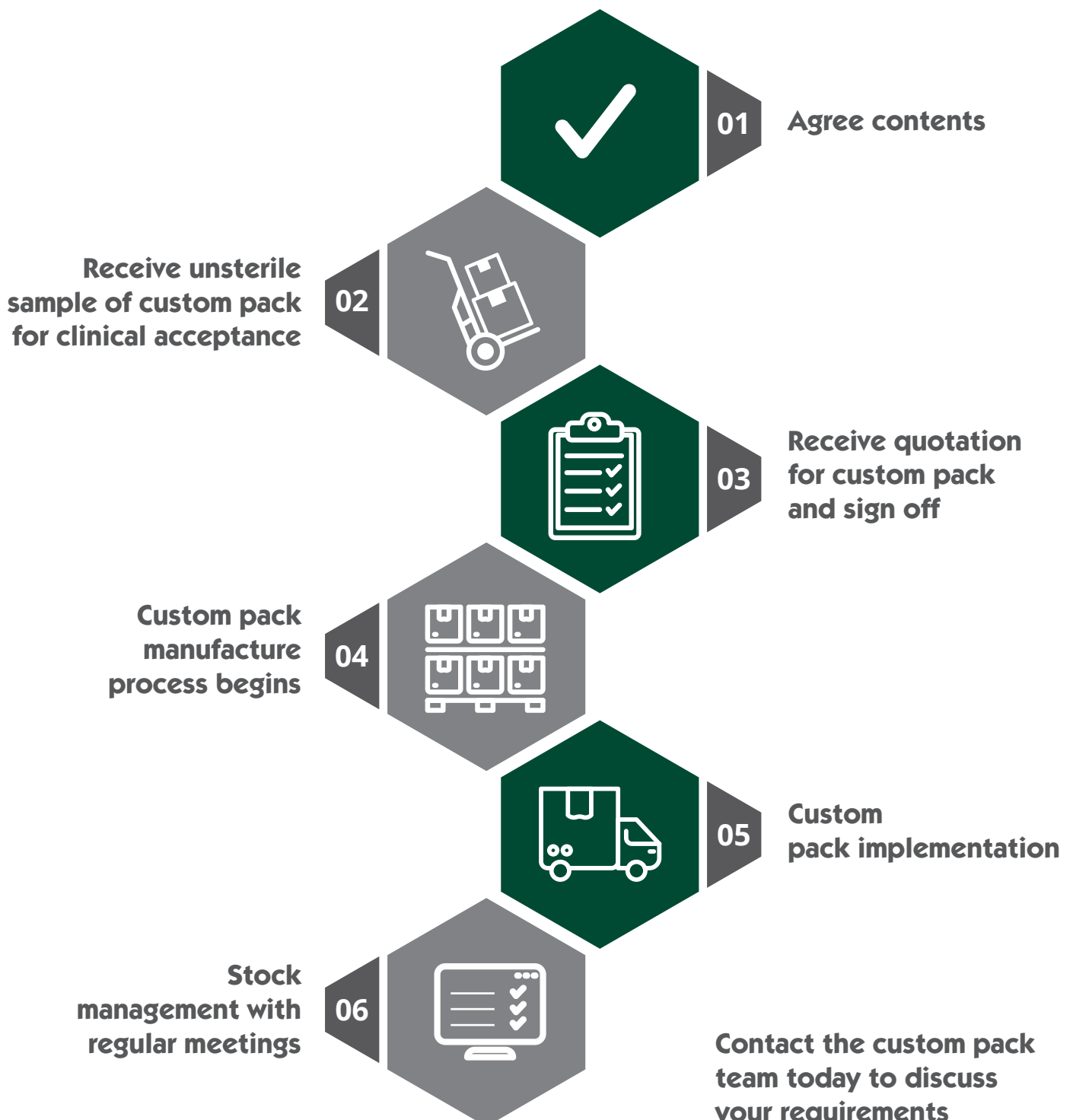

Contact the custom pack team today to discuss your requirements
01793 746904
veterinary-UK@vygon.com

“The custom packs are designed to our own specifications which is great for us.

It saves us so much time all round and is no more expensive than buying everything separately. The boxes arrive sterilised and ready to use. We get them in as and when we need them and it's just a great service.

As well as saving us time in ordering the products separately from different suppliers, we also save wear and tear on the autoclaves. We want to spend as much of our time with our patients and their owners as we can so the packs are just perfect for what we need.”

Nicola Casement

Referral Veterinary Nurse, Manchester Veterinary Specialists

Why choose custom packs procedure packs?

From concept to implementation we will work with you to develop a high quality pack to meet your individual needs, giving you peace of mind throughout:

- **Standard practice**
- **Sterilisation**
- **Warehousing & distribution**
- **Implementation & training**
- **Stock management**



Reduce waste



Save money



**Reduce
carbon footprint**



**Save on
set up time**



ONCOLOGY

Qimono Qimooctopus

Qimooctopus with Qimo female and Qimo male is a multi-purpose catheter accessory, a double lumen IV extension tube with needleless connectors allowing intermittent injections or continuous infusions of fluids or medications.

Code	Box qty
007231221	1



Qimono Qimocap

Qimocap is a vial access cap fitted with Qimo female.

Code	Size (mm)	Box qty
00721901	14	50
00721902	20	25



Qimono Qimospike V

Qimospike V is a vented spike for vials fitted with Qimo female.

Code	Box qty
007217040	25



Qimono Qimomale

Qimomale is a multi-purpose, closed, needle-free, IV connection device which has been designed for the preparation and safe administration of hazardous drugs. Qimomale has a very easily disinfectable membrane (disinfection must be carried out before and after use) which automatically closes when Qimo female is disconnected.

Code	Box qty
00721091	50



Qimono Qimofemale

Qimofemale is a multi-purpose, closed, needle-free, I.V. connection device.

Code	Box qty
00721002	50



ONCOLOGY PPE

Non-sterile chemotherapy gowns

Code	Description	Box qty
VHU69602	Non-sterile thumb loop Chemo gown L	100 (bagged in 10s)
VHU69603	Non-sterile thumb loop Chemo gown XL	100 (bagged in 10s)

Sterile chemotherapy gowns

Code	Description	Box qty
VHU44690	Aero chrome gown XL	28

Chemotherapy gloves

Code	Description	Box qty
VHU50600M	Purple Nitrile-xtra glove XS	50
VHU50601M	Purple Nitrile-xtra glove S	50
VHU50602M	Purple Nitrile-xtra glove M	50
VHU50603M	Purple Nitrile-xtra glove L	50



NUTRISAFE2® - SAFE ENTERAL FEEDING SYSTEM

nutrisafe2® is Vygon's complete and unique system for secure enteral delivery in small animals.



Safe

Nutrisafe2 is incompatible with standard 6% luer taper devices. The lockable non-luer system reduces the risk of misconnection errors with no change to normal practice. The Nutrisafe2 connector and the plastic cap tether have been designed to ensure it is safe for even your smallest patients.



Secure

The lockable system provides additional security. Accidental disconnection is virtually impossible. To disconnect the system, you do not need to use force, you simply unscrew the device.



Confident

Colour-coded for easy recognition and size identification. The Nutrisafe2 range features purple hubs and size related colour-coded caps, reducing the risk of accidentally attempting to connect devices intended for use with IV lines. The unique Nutrisafe2 secure lockable non-luer system gives you peace of mind that your patients are safe from misconnections and wrong route drug errors.



Comfortable

The smooth shape was chosen specifically to eliminate sharp edges that could cause discomfort and markings on delicate skin.



A Complete Solution

Vygon offers a wide range of syringes, extension tubes, drawing-up straws, and accessories. Numbered graduations every centimetre aid placement and monitoring, with a rounded tip providing atraumatic insertion.

With continual support offering team training and technical assistance to support best practice within your practise.

ENTERAL FEEDING

nutrisafe²® PVC nasogastric feeding tubes

DEHP-free PVC Feeding Tube – Short-term use			
Code	Size (Fr)	Length (cm)	Box qty
000364042	4	40	25
000361042	4	50	50
000363042	4	75	25
000364052	5	40	25
000361052	5	50	50
000363052	5	75	25
000364062	6	40	25
000361062	6	50	25
000363062	6	75	25
000362062	6	125	50
000364082	8	40	25
000361082	8	50	25
000363082	8	75	25
000362082	8	125	25
000363102	10	75	25
000362102	10	125	25
000363122	12	75	25
000362122	12	125	25

nutrisafe² drawing-up straws

Drawing-up straws allow the filling of nutrisafe² syringes with feed or medication from all sizes of containers and bottles.


Code	Length (cm)	ID (mm)	OD (mm)	Box qty
000817152	15	2.0	3.3	50
000817153	15	3.0	4.2	50



nutrisafe² PUR nasogastric feeding tubes


The radiopaque feeding tubes come with size colour-coded cap and are suitable for medium-term use.

PUR Feeding Tube - Medium-term use			
Code	Size (Fr)	Length (cm)	Box qty
001361042	4	50	50
001363042	4	75	25
001361052	5	50	50
001363052	5	75	25
001362052	5	125	25
001361062	6	50	50
001363062	6	75	25
001362062	6	125	25
001361082	8	50	50
001363082	8	75	25
001362082	8	125	25
001361102	10	50	50
001363102	10	75	25
001362102	10	125	25
001363122	12	75	25
001362122	12	125	25
001362142	14	125	25
001362182	18	125	25

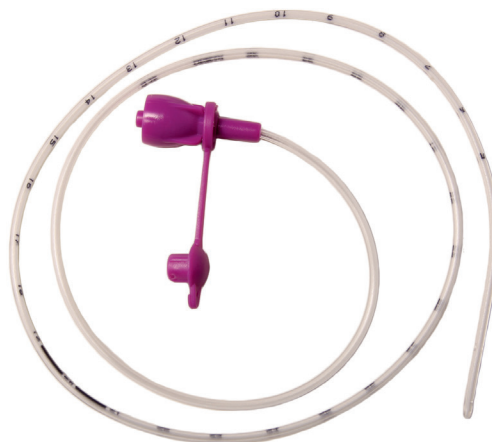

nutrisafe² PUR feeding tubes with Stylet.

PUR Feeding Tube - Medium-term use			
Code	Size (Fr)	Length (cm)	Box qty
0013631421	14	75	25
0013631621	16	75	25
0013620822	8	125	25

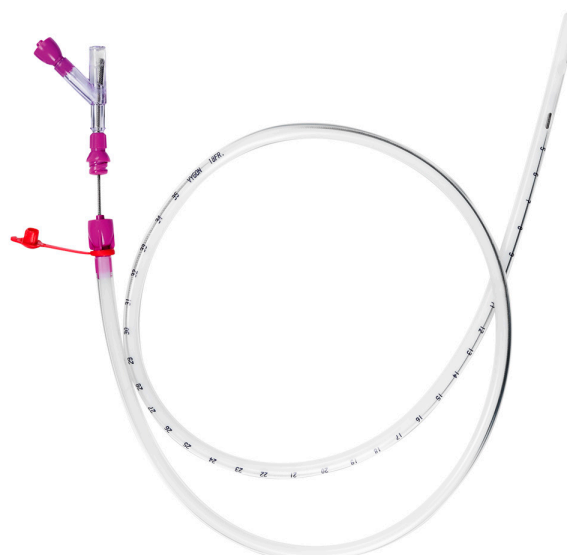

nutrisafe2 extension tubes

This extension tube gives extended connections to pumps. It can only be used with feeding tubes featuring a nutrisafe2 connection.

Code	Length (cm)	Length (mm)	Box qty
000368102	100	1.5-2.5	50
000817153	150	3.0-4.00	50



nutrisafe2 syringes

nutrisafe2 syringes have a purple plunger to aid instant identification as enteral devices. For safety they only connect with the other items in the nutrisafe2 range. They are graduated and available in a wide range of sizes.

Code	Description	Box qty
001015014	nutrisafe2 1mL syringe	100
001015024	nutrisafe2 2.5mL syringe	100
001015054	nutrisafe2 5mL syringe	100
001015104	nutrisafe2 10mL syringe	100
001015214	nutrisafe2 20mL syringe	80
001015604	nutrisafe2 60mL syringe	60



THEATRE

Mini Red-O-Pack

120mL container in PE/PVC for drainage of small volumes after surgery. Supplied with Redon-type suction drain, Alene needle, and tape.

Code	Description	Suction drain size (Fr)	Container (mL)	Box qty
VKRAP0614101251	Mini Red-O-Pack	6	120	1
VKRAP0814101251	Mini Red-O-Pack	8	120	1



Drape towels with adhesive

Two-layer general purpose sterile adhesive-edge drapes available in a variety of sizes. Adhesive edges provide an effective seal, preventing contamination of the surgical field from beneath the drape.

Code	Description	Size (cm)	Box qty
00412603	Drape towels with adhesive	80x90	50
00412605	Drape towels with adhesive	90x100	50
00412905	Drape towels with adhesive	180x200	20



Universal drape set

Set of four sterile adhesive-edged two-layer drapes suitable for many veterinary procedures including laparotomy. Ideal for small to medium sized veterinary patients. Suitable drape set for ventral midline ovariohysterectomy in the bitch.

Code	Description	Size (cm)	Box qty
00463912	Adhesive-edged utility drapes x4	70x75	30
004029023	Drape pack with trolley drape: x1 trolley drape – 140x190cm x2 Drape Towel with adhesive - 90x100cm x2 Drape Towel with adhesive – 30x40cm	-	18

Gowns

Sterile Surgical Gowns made from SMS material with hand towels and raglan sleeves, book folded for aseptic donning of gown.

Code	Description	Box qty
VHU11783	Sterile Surgical Gown Large with hand towels	40
VHU11784	Sterile Surgical Gown Extra Large with hand towels	40
VHU90002	Sterile Surgical Gown Small with hand towels	36
VHU90112	Sterile Reinforced surgical gown Large with hand towels	32
VHU90142	Sterile Reinforced Surgical gown extra large with hand towels	30
VHU90143	Sterile Reinforced Surgical gown Extra Large – Extra Long with hand towels	32



Fenestrated drape towels

Code	Description	Size (cm)	Box qty
00344545	Fenestrated towel with adhesive	45x45	80
00347575	Fenestrated towel with adhesive	75x75	30
00347590	Fenestrated towel with adhesive	75x90	30
003410100	Fenestrated towel with adhesive	100x100	30



Trolley / table drapes

Code	Description	Size (cm)	Box qty
00412822	Trolley/Table cover with adhesive strip	120x140	65
00412825	Trolley/Table cover with adhesive strip	140x190	42
00412826	Trolley/Table cover with adhesive strip	160x240	33



Surgical scrub brushes

Code	Description	Box qty
00090200	Chlorhexidine surgical scrub brush	25
00090100	Iodine surgical scrub brush	25
07090300	Dry surgical scrub brush without nail pick	25



Bulb syringe

Code	Description	Size (mL)	Box qty
00045100	Bulb syringe	100	30



Thoracic trocar and drain

Transparent radiopaque PVC thoracic trocar and drain with 2cm graduations and two lateral eyes. Handle ensures controlled insertion.

Code	Description	Size (Fr)	Length (cm)	Box qty
00062612	Thoracic trocar and drain	12	25	20
00062616	Thoracic trocar and drain	16	25	20
00062620	Thoracic trocar and drain	20	40	20
00062624	Thoracic trocar and drain	24	40	20
00062628	Thoracic trocar and drain	28	40	20



TIVA sets

Total intravenous anaesthesia from pump to patient.

Vygon codes	Number of lumen	Detachable extension line	Gravity line		High flow drug line		Low flow drug line		Box qty
			Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	
VAN7300A	Triple	No	8.0	0.5	10	0.6	12	0.3	50
VAN7200A	Double	No	8.0	0.5	10	0.6	-	-	50

- Longer options are available



Features & benefits



Anti-syphon valve

One-way valves that require a far greater driving pressure to open, minimising the risk of downloading.




Anti-reflux valve

Minimises the risk of back-tracking and enhances patient safety.



Halkey Roberts clamp

Secure, easy-to-use clamps, which enable one-handed operation during medication changes.

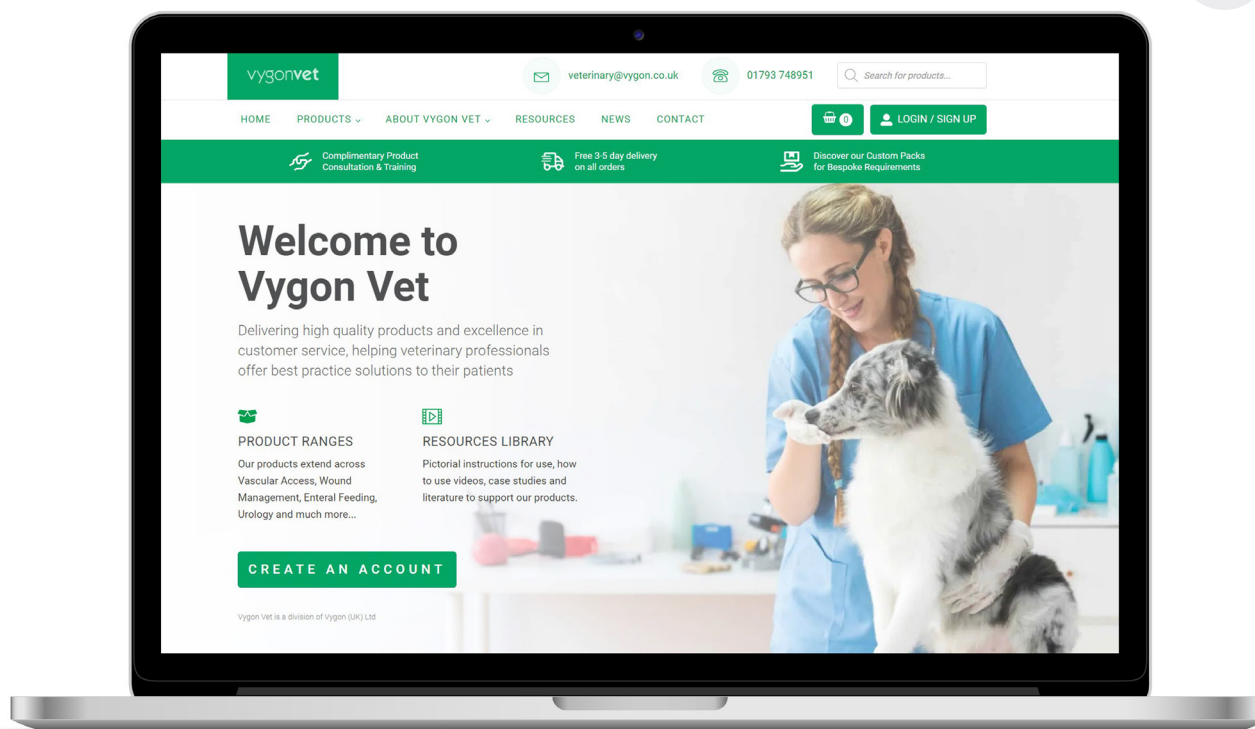

Have you registered?

Sign up for an online account today and reap the benefits



SAVE TIME

The site features all our products for instant online purchases.



ACCOUNT SPECIFIC PRICING

Having an account with a Vygon Vet can save you up to 20% on certain products.



REPEAT ORDERS WITH EASE

It takes just one click to repeat orders from previous purchases.

Take a look  vygonvet.co.uk

OUR DISTRIBUTORS

You can also buy through our distributor partners:



This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

For further information, please contact: info@vygon.com

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature. This brochure has been printed on responsibly sourced and sustainable material. To help us reduce our carbon footprint all of our literature is available electronically either from your Product Specialist, on our website or by emailing the address above.

Vygon (UK) Ltd, The Pierre Simonet Building, V Park, Gateway North, Latham Road, Swindon, Wiltshire SN25 4DL
Reception: +44 (0)1793 748800 www.vygon.co.uk

  [@vygonuk](https://www.youtube.com/vygonuk)

